Priming **Embra® Gravity set,** vented with semirigid container or bag


REF: 10400.

Priming recommendation to prime Embra® Gravity set when using semirigid container or infusion bag.

Please note! See the package for complete IFU instructions for use, this is only a guide for users. DO NOT OVER-TIGHTEN connections. DO NOT USE any instruments to tighten connections.Tip: Glass and non fully collapsible containers: Keep the vent closed during spiking and priming. Only open the vent flap when venting is required.


1

Unwrap the IV-set and check whether the air vent cap is closed and the roller clamp is open.

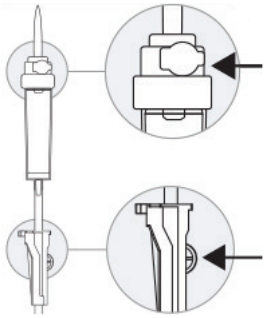


2

Change position of the roller clamp and control that the tubing is unaffected and close it.

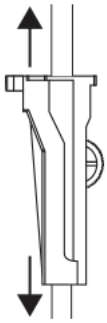


3

Spike the semirigid container/ bag. Invert the container/bag and hang it from the IV pole or similar.

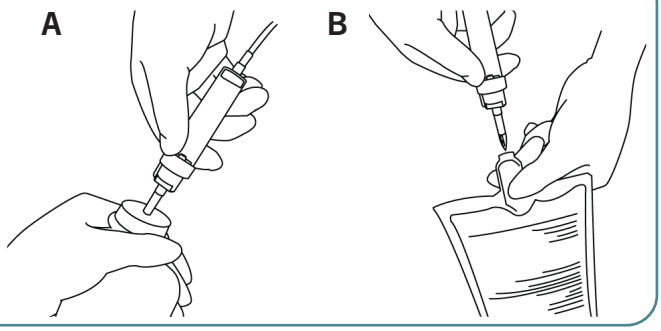


4

Note!
Attachment of the tube to the holder before squeezing the drip chamber.

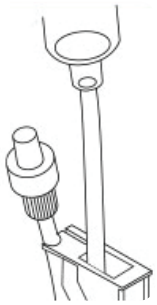

5

Squeeze gently the drip chamber 3-5 times to fill half of the drip chamber. (Too aggressive squeeze can lead to mixed air fluid below the drip chamber.)

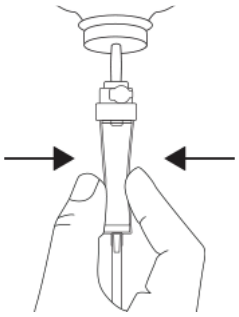

Evercare Medical – developing medical supply products for a positive difference in the European healthcare sector.

info@evercare-medical.com
www.evercare-medical.com

evercare® medical | Part of Asker Healthcare Group

Priming **Embra® Gravity set,** vented with semirigid container or bag


REF: 10400.

Priming recommendation to prime Embra® Gravity set when using semirigid container or infusion bag.

Please note! See the package for complete IFU instructions for use, this is only a guide for users. DO NOT OVER-TIGHTEN connections. DO NOT USE any instruments to tighten connections.Tip: Glass and non fully collapsible containers: Keep the vent closed during spiking and priming. Only open the vent flap when venting is required.


6

Open the roller clamp until the complete tubing is primed.

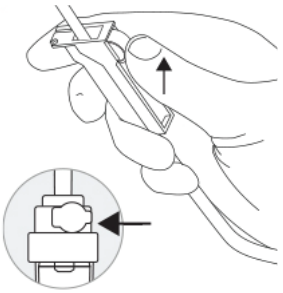

7

Close the roller clamp.

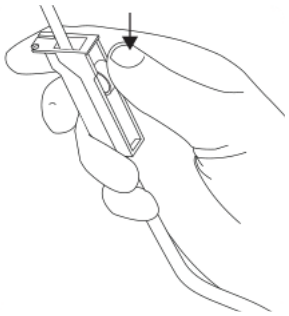

Evercare Medical – developing medical supply products for a positive difference in the European healthcare sector.

info@evercare-medical.com
www.evercare-medical.com

evercare® medical | Part of Asker Healthcare Group